

## Sound proof cabinets for ultrasonic baths

Sonation USBB series



**-30dB(A)  
Over 90% reduction in noise**

Easy access


Space-saving

Robust construction

Active ventilation



- Noise reduction up to 30dB(A) which is a drop in noise by over 90%. The cavitation noise is practically eliminated.
- Available in four different sizes for mostly all ultrasonic baths.
- The cover can be completely opened, so the removal of the basket is very easy to handle.
- There is no additional space needed to open the cover, so the box can be placed and operated under the extractor hood or under lab cabinets.
- Control panel stays visible. Thanks to the transparent front, the progress of the washing program can be observed anytime.
- Equipped with fast access openings on the front and both sides to operate the ultrasonic bath or to drain the cleaning agent.
- Integrated fan to remove vapors.
- Oil and water repellent and flame-retardant acoustic foam.
- Powder coated, robust and easy cleanable surface.
- Cushioned opening mechanism for smooth and easy handling.

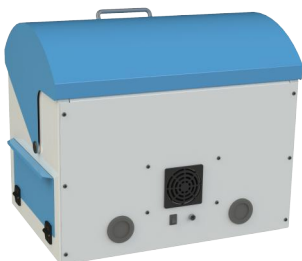

### Contact:

Sonation GmbH  
Alte Schulstr. 39  
88400 Biberach  
Tel.: +49 (0) 7351 30038-39  
E-Mail: mail@sonation.com

[www.sonation.com](http://www.sonation.com)

## Technical data:

	<b>USBB-S</b>	<b>USBB-M</b>	<b>USBB-L</b>	<b>USBB-XL</b>
Outer dims. (WxDxH)[mm]:	512x350x425	557x410x545	607x460x633	782x610x671
Inner dims. (WxDxH)[mm]:	454x269x316	479x329x415	529x379x496	704x529x516
Weight [kg]:	15	20	25	30

## Compatibility list:

The following list should help you to choose the right model. Please note, however, that the dimensions of ultrasonic baths by the manufacturers are often specified without drain taps, handles or other attachments. It is therefore essential to measure whether the corresponding model actually fits into the sound insulation housing.

### Advantage-Lab

**USBB-S:** AL04-01, AL04-02, AL04-03, AL04-04, AL04-06, AL04-12

**USBB-XL:** AL04-30-230, AL04-45, AL04-90

### BANDELIN

**USBB-S:** DT31(H), DT52(H), DT100(H), DT102H(H-RC), DT103H, DT106, DT255(H)(H-RC), DT510(H)(H-RC), DT512H(H-RC), RK31(H), RK52(H), RK100(H)(SH), RK102H(CH), RK103H(CH), RK255(H), RK510(H), RK512H, RK255CH, DK102P, DK255P, DK512P

**USBB-M:** DT514, DT514H, DT514BH, DT514BH-RC, DT106, RK513, RK514, RK514H, RK514BH, RK512CH

**USBB-L:** Bactosonic 14R, Sonorex Digital 10P, DT106, DT156, DT156H, DT1028, DT1028H, RK106, RK156, RK156BH, RK255H, RK1028, RK1028H, RM16

**USBB-XL:** RM40, RM40-U, RM40-H, RM40-UH

### BRANSON

**USBB-S:** Branson 200, 1510, 1800, 2510, 2800, 3510

**USBB-M:** Branson 3800

**USBB-L:** Branson 5510, 5800

**USBB-XL:** Branson 8510, Branson 8800

### Elma

**USBB-S:** Elmasonic S10/H, S15/H, S30/H, S40/H, S60/H, S90/H, S100/H, S120/H, S130/H, P30H, P60H, P120/H, Denta Pro, xtra TT 30H, xtra TT 60H, S50R, P120H, TI-H-5, T710DH, T760DH, T490DH

**USBB-M:** Elmasonic S180/H, P180H, xtra TT 120H, xtra TT 150H

**USBB-L:** Elmasonic xtra TT 200H, TI-H-10, TI-H-15, T680DH

**USBB-XL:** Elmasonic S70/H, S80/H, S150/H, S300/H, S450/H, S900/H, P70H, P300H, TI-H-20

### Grant Instruments

**USBB-S:** XUB5, XUB12, XUBA1, XUBA3

**USBB-M:** XUB18

**USBB-L:** XUB10, XUB25

### SOLTEC

**USBB-S:** Sonica 1200, 2200, 2400, 3200, 3300

**USBB-M:** Sonica 4200, 4300

**USBB-L:** Sonica 3200 L

**USBB-XL:** Sonica 5200, 5300, 45L EP, 90L EP, ATC-EP

### VWR

**USBB-S:** USC 100, 200, 300, 500, 600, 1200

**USBB-M:** USC 900, 1700

**USBB-L:** USC 2100, 2600

### US-Cleaner

**USBB-L:** 2800, 4500

### Contact:

Sonation GmbH  
 Alte Schulstr. 39  
 88400 Biberach  
 Tel.: +49 (0) 7351 30038-39  
 E-Mail: mail@sonation.com